

INSTRUCTIONS FOR USE

Congratulations on purchasing your PHALLOSAN® forte! PHALLOSAN® forte has been being used around the world with much success for 18 years. Prescription PHALLOSAN® forte is reimbursed by the health insurance companies. Our team consisting of urologists and medical technicians developed what we think is the most effective and innovative product of its kind, PHALLOSAN® forte.

Indications:

- PHALLOSAN® forte is recommended as a treatment for
- » penis enlargement-penis extension
 - » congenital/acquired penile curvature and induratio penis plastica (IPP)
 - » erectile dysfunctions
 - » post-surgical follow-up treatment after prostate surgery
 - » training the cavernous body of the penis at any age

Package contents:

1. Elastic belt, which is strapped around the waist incl. foam-rubber ring (A), into which the penis is inserted all the way to the base. The buckle (B), with which you can adjust the length of the belt individually and adapter (C) on the belt, which is connected to the suction bell.
 2. a) Suction bell in size S with PHALLOSAN® forte – S-sleeve condom
 2. b) Suction bell in size M with PHALLOSAN® forte – M-sleeve condom
 2. c) Suction bell in size L with PHALLOSAN® forte – L-sleeve condom and assembled tension clip (3)
 3. Tension clip with tension spring attachment to adjust the intensity of the tension using a traffic light colour code system with adapter for the belt.
 - 4./5. Protector cap to protect glans against strong tension settings.
 6. Suction ball with 3-way valve
 7. Measuring template
- Instruction Guide
Traveller Bag with 4 extra sleeve condoms (1xS, 2xM, 1xL)



Preparation

1. Belt



Spread the orthopaedic PHALLOSAN® forte belt out in front of you and adjust it to its longest length using the buckle (B).

2. Selection of the suction bell/sleeve

Determine your correct suction bell size using the enclosed measuring template (7). Pull the flaccid penis away from your body and take measurement on the penis shaft.



The bell in size M (2c) is already connected to the sleeve size M (2c) and the tension clip with tension spring attachment in the set.

Should you require a different size, please remove the tension clip (3).

Now, please take your appropriate bell and adapt it using the suction ball; it only snaps into the position with the arrow symbol (↕). The appropriate sleeve (2) is slipped over the bead of the bell and the tension clip (3) is reattached.



When attaching to the bell, it must be observed that the bead of the bell snaps into the groove on the sleeve on the thick end of the sleeve. If the sleeve is pulled too far over the bead of the bell, the sleeve can rupture.





Top: Wrong, pulled too far over the bell and not powdered!

Bottom: Right and powdered.

Please observe that the sleeve loses its coating after several days. Therefore, you should purchase baby powder and apply this to the outer surface of the sleeve. It forms a protective coat, which protects the sleeve against friction, adhesion and tearing.

Clamp it the suction bell between your knees. With a slight pull, carefully roll the sleeve condom 0.40 inch to the edge of the suction bell and pull it behind the bead of the suction bell.

Please be careful and do not damage the sleeve with sharp fingernails or the like. Do not pull hard when rolling back.

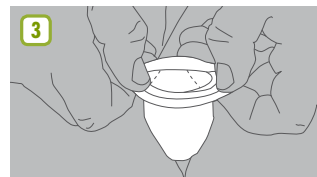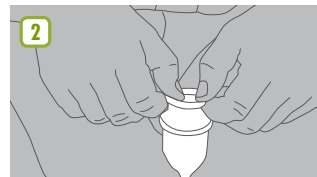





Special features of the new PHALLOSAN® forte sleeves

The PHALLOSAN® forte sleeves consist of medical silicone, thus ruling out the risk of irritation or allergic reactions. This innovation not only offers you more comfort but also provides a uniquely snug fit around the penis. A part of it is approx. 3 mm thick, that is where it is attached to the bell; its thinnest area is 0.3 mm thick, to smoothly adapt to the shaft of your penis.

Your perfect suction bell is now ready for use!

3. Protector cap

Please remove the protector cap (4/5), which fits any penis size. Pull the protector cap apart using both hands and slip it over the glans. It should fit snugly.

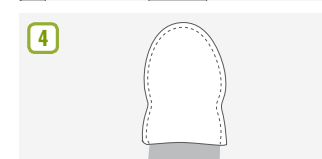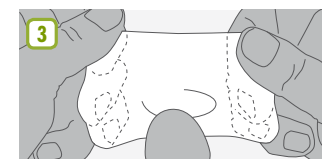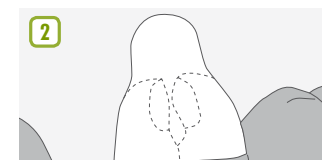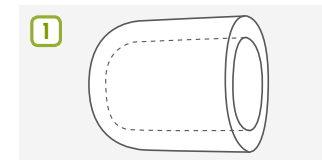

After fitting the cap, it is important to make sure that as little air as possible remains at the tip of the glans. It can be helpful to slip the protector cap over the penis when the latter is in a slightly erected state. The protector cap should protect the glans against the effects of possible excess tension.

Experienced users can usually do away with the protector cap after a few weeks.

All necessary preparations for the application have now been completed.

Application:

1. Now, insert the penis including the protector cap into the suction bell and roll the sleeve back from the bead of the bell over the penis.



2. Turn the green protruding marker on the valve of the suction ball (6) to the suction symbol and press the suction ball several times with your fingers. By doing so, the glans is drawn further into the bell. The sleeve smoothly adapts to the shaft of the penis.



3. Now, turn the green protruding marker on the valve of the suction ball (6) to the lock symbol, to maintain the micro-vacuum.



4. Now, please take the belt previously spread out before you and slip the white foam-rubber ring (A) over the complete suction bell apparatus until it rests against the root of the penis.

5. Now, pull the end of the belt around your waist.

6. Turn the open end of the tension clip (3) towards your body and connect both eye-lets (C).



7. Tighten the fit of the belt using the belt buckle (B) until you can feel that your penis is being stretched and the traffic light colour code system indicates green. Please make sure that the belt slants upwards in the groin area. To do this, pull the belt upwards around the buttocks towards the waist. If the belt slips to the fold of the buttocks, i.e. parallel to the penis, the tension will be further increase when sitting, which may be unpleasant.



8. You can now carry on with your everyday business or also wear PHALLOSAN® forte while sleeping, given that your blood circulation is not restricted at any time.

Ending the application

First, please take the belt off and turn the green valve of the suction ball to the arrow symbol. The vacuum escapes. Please roll the sleeve carefully back towards the hard plastic cup and remove the complete suction device. **Do not stretch the sleeve with your fingers, to remove the vacuum, as it will otherwise rip immediately!** It will be easily removed when simply rolling the sleeve upwards.

Special technical features

Put simply, this apparatus is a suction bell with a seal. Therefore, it is not important how far the sleeve reaches over the penis shaft. The sleeve is not designed to transmit force to the penis. This solely takes place by the micro-vacuum formed in the bell.

Should you feel that the vacuum has slightly decreased while wearing the apparatus and the glans has possibly slipped too far out of the bell, this can be the result of several factors. In any case, air is penetrating. This can take place via the penis shaft if pubic hair was not diligently removed or a sleeve, which is too wide, was selected.

As far as necessary, please remove possible pubic hair from the penis. This can be done by shaving or using a mild hair removal cream.

Remedy

Turn the green protruding marker on the valve of the suction ball away from the lock symbol to the suction symbol. In this manner, draw the glans back into the bell and close the valve by turning the green protruding marker back towards the lock symbol.

If the glans is still slipping out of the bell, the cause is definitely the sleeve, which is apparently too wide.

You have the option of adjusting this individually to your size and in doing so, achieving optimum efficiency. For example, if you feel that the M bell is the correct size but the sleeve is too large for the penis shaft, you can place the S sleeve on the M bell. You can do the same with the L bell and the M sleeve.

Under no circumstances can you place the S sleeve on the L bell. This can cause the material to rupture.

Instructions for care:

The sleeve can and should be detached from the bell and should be washed in lukewarm water for hygienic reasons.

Keep the sleeve rolled up until you need it again.

Please never forget to apply baby powder to the outside of the sleeve after use when rolled out.

The suction bell should also be cleaned daily as described above. Do not forget to remove the suction ball first by moving the green protruding marker on the valve of the suction ball towards the double arrow. Carefully pull the suction ball out of the suction bell. The suction ball may not come in contact with water or oil because the sensitive valves would then no longer function efficiently.

Step-by-step plan for the precise adjustment of the tension of the belt to the connective tissue of the penis.



Preamble:

The body's connective tissue depends on the respective hereditary disposition. Therefore, we developed a traffic light colour code system, to enable you to adjust the belt tension precisely to your type of connective tissue. In any case, please make sure that the tension is in the lower green range in the beginning. This will help you achieve optimum results for penis enlargement in length and width without any negative side-effects.

Day 1

Only set the tension in the lower green range and do not wear PHALLOSAN® for more than 1 hour. Remove the belt and turn the green protruding marker on the valve of the suction ball towards the arrow symbol. The vacuum escapes.

Carefully roll the sleeve back towards the glans and remove the suction bell apparatus.

Now, examine the point on the penis between the opening of the urethra and the frenulum on the underside of the glans. This is the most sensitive tissue on the penis. If it is reddened or swollen, the belt tension was too tight. Suspend the application for today. If you detect changes on the penis, do not begin application again until these have subsided.

Day 2 to 8

If you were able to wear PHALLOSAN® forte on the previous day without suffering from discoloration or swelling on the points in question, you can wear PHALLOSAN® one hour longer.

Extend the period of wear every day by one hour. After 8 days you will have adjusted to wearing the apparatus.

After the adjustment period PHALLOSAN® can also be worn at night. You will notice that you can feel nubs on both symbols for pumping and locking, which enables pumping or locking without having visual contact. The sleeve expands during night-time erections; circulation is not impeded, it is improved.

You will wear PHALLOSAN® forte with enthusiasm. Your penis can become measurably longer and thicker in only a few weeks. You will experience positive side effects such as a stronger sex drive and greater pleasure.

Period of application:

PHALLOSAN® should be worn 8–10 hours per day, 6 days a week, and at any time even longer!

Side effects/Precautions

PHALLOSAN® forte is used at your own risk however, it does not pose any risks if used according to the instructions. Particularly during the first weeks of application, excess tension on the penis should be avoided.

If swelling and/or blisters occur on the foreskin or on the glans, the belt will have been set too tight. It is imperative to discontinue use until the swelling or blisters have disappeared. Then set the tension somewhat looser. Please be patient, the success of the penis enlargement process can be compared to that of an orthodontic treatment. The continuity of elongation and not the amount of tension alone is a guarantee for successful results!

Occasionally you can feel pressure pain on the penis shaft; this can result from the bulging end of the sleeve resting on the silicone rubber-foam ring and therefore, cause pain.

Remedy:

Slightly roll the end of the sleeve back from the silicone rubber-foam ring; pressure pain will no longer occur.

You can also distribute the tension of the belt over the body, to reduce tension on the root of the penis. Take hold of the belt where the suction bell is connected to it and lead the tension around the body to relieve the root of the penis. It sounds complicated but it is not; the tension is now distributed throughout the entire area of the belt.

You should not wear PHALLOSAN® forte under the following circumstances:

After consuming alcohol, after taking drugs as well as painkillers, sleeping pills and tranquilizers because your natural perception of pain may be compromised and you may not be able to promptly react to pain. Avoid using PHALLOSAN® forte in the event of infections or open wounds on the penis.

How do I measure the exact increase in the length of my penis?

Please use a foot rule. A measuring tape is not suitable. Hold your penis in a flaccid state away

from your body (90 degrees). Position the foot rule on the side of your penis (as seen looking down). Press the rule against your pubic bone.

Important

Before beginning with the application, determine your exact body weight without clothing. In the event of a possible weight gain, a layer of fat forms – in particular – in the pubic area of the pubic hair around the penis. This layer of fat reduces the distance between the tip of the penis and the pubic bone, thus falsifying (shortening) the result. Physicians call this an abdominal apron.

Repeat orders

All parts for PHALLOSAN® forte can be reordered individually from us. Should you have any other questions or suggestions, feel free to contact our experts at any time free of charge. You merely have to pay the normal telephone charges. At all times, we reserve the right to make technical modifications, which we will also make available on the Internet at

www.phallosan.com
www.phallosan.co.uk

Hotline:

Hotline +423 235 84 77
(Liechtenstein, Mo-Fr 14-16 Uhr MEZ)
service@swissana.com
for more Information:
www.phallosan.info

Your PHALLOSAN® forte Team

Read the complete instruction for use before using the product!

WARNING: The following precautions must be observed

The device **must not** be worn during bowel movements, during sports activities, during intercourse or other strenuous physical activity.

Neither should it be worn after consumption of alcohol or ingestion of pain-killers or euphorants or if there are

- any wounds or infected areas on the penis or where PHALLOSAN is attached at the body.
- If pain, numbness or discoloration of the penis head,
- or any other untoward reaction occurs the device must be loosened and removed immediately.

In case of wounds on the penis the device **must not** be worn until these have healed entirely.

In the case of surgery, **do not** use PHALLOSAN forte until your doctor tells you to.

Do not use PHALLOSAN forte if you have problems with venous problems or connective tissue weaknesses.

Do not use the device when you have any kind of piercing at the genital area.

The device **should not** be worn while moving in areas where there is an increased risk of falling and caution must be exercised when climbing stairs, walking on slippery floors etc.

The device **should not** be worn neither when handling any kind of machinery, motorized vehicles or hot objects nor when bicycle or cooking.

Do not use any creams, Vaseline and other moisture substances while wearing/using the device.

Use the device **only as we describe** in the manual, do not use any other materials than the ones we provide you or use other wearing techniques than the ones we tell you in the manual.

We do not take any responsibility and we are not liable for any damage caused through use of products or services purchased through this website, be it indirect, special, incidental or consequential damages (including but not limited to damages for loss of business, loss of profits, interruption or the like). If you have any questions regarding the terms of use outlined here, please do not hesitate to contact us at any time.

PHALLOSAN[®] forte



**Danke, dass Sie sich für PHALLOSAN® forte entschieden haben /
Thank you for purchasing PHALLOSAN® forte!**


**Holen Sie sich die PHALLOSAN® forte App gratis /
Get PHALLOSAN® forte App for free:**



Scannen Sie den QR-Code mit Ihrem Smartphone oder Tablet um die App herunterzuladen.

Scan the QR code with your smartphone or tablet to download the app.

Wir sind auch auf / Follow us on:

 Instagram — Dr.Phallosan

 Facebook — Phallosan



**Im Teamwork funktioniert es einfach besser:
Jetzt NEU der PHALLOSAN plus+.**

It just works better in teamwork: now the NEW PHALLOSAN plus+.



**Get your personal
Training Plan for
smartphone and tablet**

